


Separation of Aromatic Amines Reliasil ODS, 5 μ m

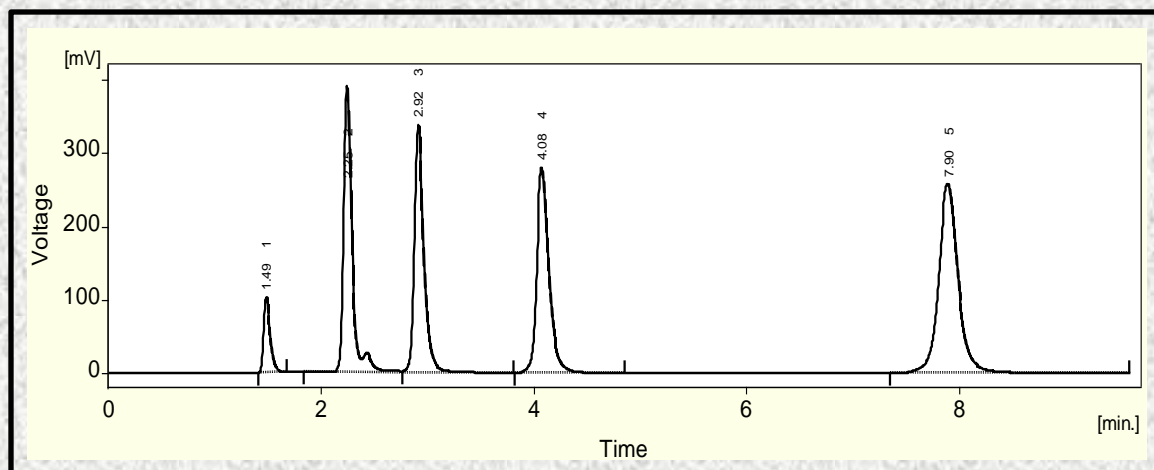

Uracil R₁	1.49 min
Aniline, R₂	2.26 min
P-Toluidine, R₃	2.92 min
4-Ethylaniline, R₄	4.08 min
N,N-Diethylaniline, R₅	7.90 min

Column	Reliasil ODS, 5 micron, 4.6 x 150 mm
Mobile Phase	60/40 Methanol/water
Flow Rate (mL/min)	1.0
Temperature (°C)	25
Detection	UV at 254 nm